



CRP FUSION • CLEARLY EXPLAINED

Rapid Drug Tests for Saliva and Urine Samples

Find the right test, use it correctly, and understand the result.

Product images: CRP 2026 catalog. Photographic scenes are illustrative.

Contents

FIND THE RIGHT TEST WITH EASE

From saliva tests and urine panels to the substance lists, the following pages guide you step by step through the range.

01 CRP Fusion Company	03
02 Guidance and Overview	04-05
03 Saliva Products	06-13

SALIVA / PRODUCTS AND USE

Alcohol Test Pen Fit2Drive	06
THC Cassette 3.5 ng/mL	07
8 + 1 Cassette	08
9-Substance Pen Industry	09
1 to 12 Substances	10
Cup 15 Substances	11
Other Products and Testing Steps	12
Understanding the Result	13



URINE / SERVICES / APPENDICES

Urine Product Range	14
Urine Panel 9 Substances	15
Contact and Selection Guide	16
Custom ALC Labels	17
Appendix A Saliva Substances	18
Appendix B Urine Substances	19-20
Price List All Product Groups	21
FAQ Products and Services	22-24

CRP Fusion - Expertise in Drug Testing and Services

BASED IN BRAȘOV SINCE 2015

CRP Fusion SRL was founded in Brașov in 2015. It is the exclusive partner of RoDiagnostics SRL, a manufacturer of saliva and urine drug tests. CRP Fusion focuses on marketing and distributing these products. In Romania, it works with Ropharma S.A. and Dasco Distribution.

PEDT - PRE-EMPLOYMENT DRUG TESTING

Since late 2025, CRP Fusion has developed this service with the PEDT department of Securitas India. According to the company, it serves the United Kingdom, Belgium, Canada, Mexico, Malaysia, and Argentina on behalf of Securitas. Expansion to Australia, South Africa, the Philippines, South Korea, and China is planned for 2026.

PARTNERS AND COOPERATION



Securitas India

PEDT Cooperation



jobselection

Job Selection

Background Screening



cFIRST

International Expansion

2027 | 50-100 CANDIDATES / MONTH

Target from early 2027

END OF 2027 | AT LEAST 200 / MONTH

Planned global expansion with cFIRST, a background screening company with offices in India, the UK, and the USA

PRODUCTS AND AUTHORIZATION

The products are certified and authorized as self-tests by ANMDMR Romania. This certification is an essential requirement for distribution in the European market.

What Do the Tests Detect?

A rapid test can screen for one substance or several substances at the same time.

SALIVA

A soft collection tip gathers saliva from the mouth. Single tests come as pens or cassettes; multiple tests come as pens, cassettes, or cups. Results usually appear within minutes, according to the product instructions.

URINE

The sample is collected in a cup. A multiple test can screen several substance groups. The exact makeup of the urine panel is not specified in the supplied catalog and must be confirmed before ordering.

IMPORTANT TO KNOW

A rapid test gives an initial indication. A reactive result should be clarified with an appropriate laboratory analysis before important decisions. A nonreactive result does not establish fitness to drive or work.

Six Ways to Choose a Saliva Test




Alcohol Pen

1 Substance



THC Cassette

1 Substance

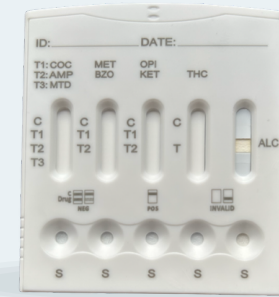

8 + 1 Cassette

8 Substances + Alcohol



9-Substance Pen

9 Substances



12-Substance Pen

up to 12 Substances



15-Substance Cup

15 Substances

The following pages explain each format using catalog information and original product images.

Alcohol Test Pen | Fit2Drive



What Does It Do?

The pen screens for alcohol in saliva. The product list includes versions with indications from 0.00 to 0.30%. These saliva-test percentages are not equivalent to breath or blood alcohol values.

How to Use It

Remove the cap, move the collection tip along the gums and tongue, and hold it under the tongue. The catalog specifies reading the alcohol indicator after 3-5 minutes. Always follow the package instructions.

PRODUCT CODE: #CRP-D062-8003-NM | also listed as a pack of three

A negative result does not certify fitness to drive. If in doubt after drinking alcohol, do not drive.

THC Cassette | 3.5 ng/mL



THC = Active Compound in Cannabis

The test detects THC in saliva. Its result does not show the degree of impairment or whether a legal blood limit is met.

PRODUCT CODE: #CRP-D022-P8002-NM

The source states “120 ml” for three drops; that plainly inconsistent volume was omitted.


What Does 3.5 ng/mL Mean?

This is the cutoff stated in the catalog for the saliva cassette. The same number used for a legal blood-serum limit refers to a different measurement in a different sample. A blood value cannot be calculated from a saliva result.

Simple Procedure

Collect saliva with the sponge, press it into the tube, and apply three drops to the cassette. Read the result within the stated 3-10-minute window. The specific product instructions take precedence.

Fit2Drive 9 | 8 + 1 Cassette

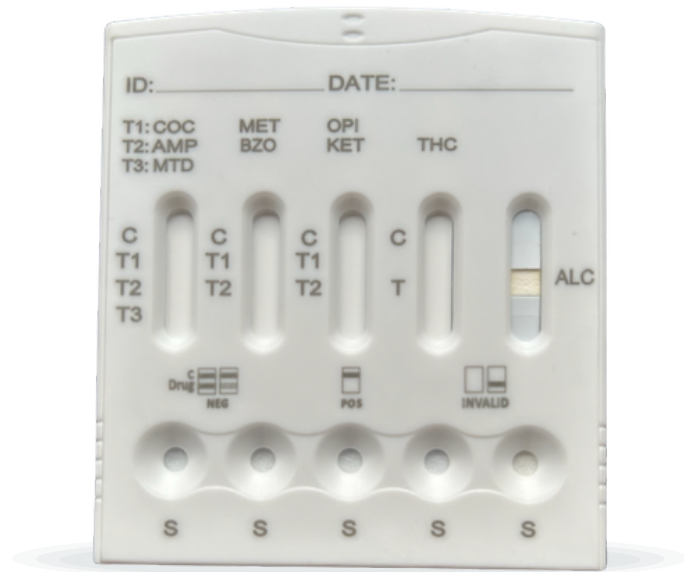

Who Is It For?

One cassette screens several substance groups in a saliva sample. The catalog presents it as a Fit2Drive version. It gives a screening result, not a fitness-to-drive assessment.

PRODUCT CODE: #CRP-D000-8095-A-NM

OP is the code used in the 2026 product list. The precise target and all cutoffs are stated in the product specification.

COC

Cocaine

AMP

Amphetamine

MTD

Methadone

OP

Opiates

MET

Methamphetamine

BZO

Benzodiazepine

KET

Ketamine


THC

Cannabis

ALC

Alcohol

9-Substance Pen | Two Workplace Versions

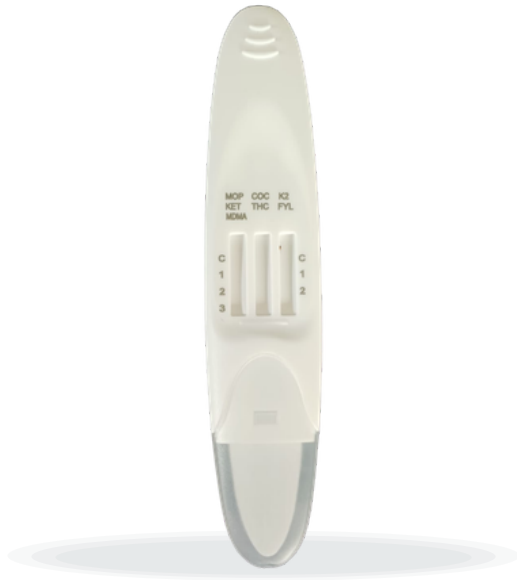

Office Version

MOP • KET • MDMA • COC • AMP • PCP • MET • THC • LSD. The pen screens nine catalog-listed targets at the same time.

Production Version

K2 • MTD • TML • MOP • BZO • THC • FYL • BAR • OXY. A different nine-target set for the applicable testing policy.

How Is the Pen Used?

Remove the cap, move the collection tip along the gums and tongue for about 60 seconds, and hold it under the tongue. Read the result according to the instructions; the catalog states approximately three minutes at the earliest.

PRODUCT CODE: Office: #CRP-D000-8093-D-OF-NM

PRODUCT CODE: Production: #CRP-D000-8093-D-PR-NM

Letter codes identify panel targets. Check the precise substance definitions and cutoffs for the specific product version.

One Substance or Up to 12 in One Test



Multiple Pen: 2-12 Substances

The catalog describes configurable saliva pens screening 2 to 12 targets. It does not specify one fixed 12-target configuration. Select the substance groups before ordering and obtain written confirmation of the panel, cutoffs, and product code.

ALSO LISTED AS SINGLE-TARGET PENS

COC

Cocaine

20 ng/ml

KET

Ketamine

50 ng/ml

MDMA

Ecstasy

50 ng/ml

THC

Cannabis

12 ng/ml

LSD

LSD

10 ng/ml

These are saliva cutoffs from the product list. The single-target THC pen (12 ng/mL) differs from the 3.5 ng/mL cassette.

Cup | 15 Substances in One Test



15 SUBSTANCES IN ONE TEST PRODUCT

MOP Opiates	MET Methamphetamine	OXY Oxycodone
COC Cocaine	AMP Amphetamine	LSD LSD
BZO Benzodiazepine	MDMA Ecstasy	FYL Fentanyl
KET Ketamine	ZOP Zopiclone	MTD Methadone
BAR Barbiturates	TML Tramadol	THC Cannabis

15 Test Strips in One Cup (CUP)

Collect saliva with the sponge and place it in the cup. Close the lid and mix as instructed. The catalog specifies reading after 3-10 minutes. Interpret each test strip separately.

PRODUCT CODE: #CRP-D000-8157-C-F2D-NM

The cup screens substance groups; it does not certify fitness to drive. A reactive result needs confirmation.

Other Products and Testing Steps

Fit2Work Office | 12 Targets

Example configuration: MOP, MET, COC, AMP, PCP, LSD, MDMA, KET, MEP, MDPV, alpha-PVP, and THC.

Fit2Work Production | 12 Targets

Example configuration: MOP, OXY, MDMA, BZO, FYL, CFYL, K2, ZOP, MTD, BAR, TML, and THC.

FIVE SIMPLE STEPS

1

Open

Check the seal and expiration date.

2

Collect

Collect saliva with the soft tip.

3

Insert

Place the collection part in the cup as instructed.

4

Wait

Observe the time window on the package.

5

Read

Check the control and test lines for each field.

Office: #CRP-D000-8127-C-F2WO-NM • Production: #CRP-D000-8127-C-F2WP-NM

Before use, check each cup configuration and its cutoffs against the current product specification.

What Does the Result Mean?

NONREACTIVE

The control is valid. The substance was not detected at this cutoff. This does not establish abstinence or fitness.



C 
T 

REACTIVE

The screening indicates a preliminary reactive result. Obtain an appropriate confirmation before any consequential decision.



C 
T

INVALID

The control line is absent or the display cannot be interpreted. Repeat the test with a new device as instructed.



C
T 

C = control line; T = test line. This diagram applies to drug strips; alcohol indicators may differ.

Why Do Cutoffs Matter?

A cutoff defines the concentration at which a test is designed to react. Saliva and urine are different sample types. Even identical numbers in saliva and blood cannot be compared directly.

For prescribed medicines, review the result and medication confidentially with a qualified professional.

Displays vary by product. This diagram does not replace the instructions for the specific test.



A SEPARATE PRODUCT CATEGORY

Rapid Drug Tests Using Urine

Multiple panels with their own collection procedure and product specifications.

The background scene is illustrative and does not show a specific CRP product version.

Urine Panel | 9 Substances

The supplied product label shows a specific CRP multiple urine test.



NINE DRUG TARGETS

One urine sample, nine test fields:

OPI

Opiates

PCP

Phencyclidine

COC

Cocaine

BAR

Barbiturates

BZO

Benzodiazepine

AMP

Amphetamine

MTD

Methadone

TCA

Tricyclic
Antidepressants

THC

Cannabis

PRODUCT CODE: CRP-D000-1094B-NM

OTHER URINE PANELS

CRP also lists 12- and 20-target versions. Request the current target list, cutoffs, reading time, and product images for these versions separately.

The 9B label comes from the supplied image. Check the current package instructions and product specification before use.

Contact and Product Selection Guide

How to Order the Right Product

Specify the sample type (saliva or urine), the desired substance groups, the number of targets, and the intended setting. Request the current instructions, cutoffs, and product code for the exact version.

CRP FUSION

Phone Romania: +40 729 642 999

Phone International: +40 748 359 283

Phone Germany: +49 174 175 4896

contact@crpfusion.ro

www.crpfusion.ro

9 Alexandru Odobescu Street, 500271 Braşov, Romania

SOURCES AND EDITORIAL NOTES

Saliva product names, codes, panels, and images: “CS RO_CATALOG CRP Fusion_2026.pdf,” pages 1–9. Urine product 9B: separately supplied product label. Other urine panels require current data sheets.

General guidance on rapid tests: FDA, “Drugs of Abuse Home Use Test”; MedlinePlus, “Drug Testing.” German THC blood-serum limits: § 24a StVG; a saliva value is not equivalent.

Rapid tests are screening tools. The manufacturer’s instructions for the specific product take precedence.

As of September 2026 • Images from the supplied catalog; photographic scenes are illustrative.

Custom Labels for the ALC Self-Test

The labels can also be used on boxes and outer packaging.




DESIGN EXAMPLES FROM DIFFERENT PROJECTS

CRP Fusion designs and produces custom labels in-house and supplies the chosen label on a plastic pouch or in a box. The customer determines the number of tests in each box.

LABEL SIZE: 154.5 × 47.5 mm




Echipa

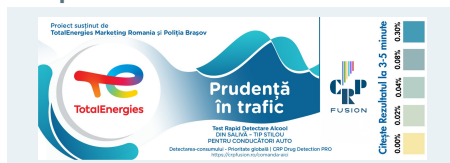


TotalEnergies




LiT

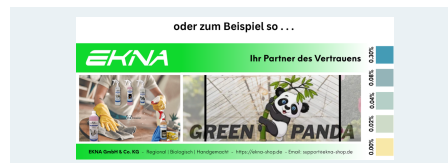


EKNA

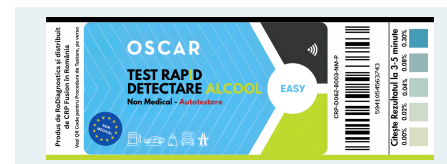

Oscar Gas

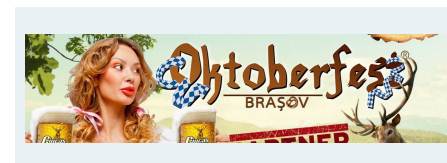

Oktoberfest

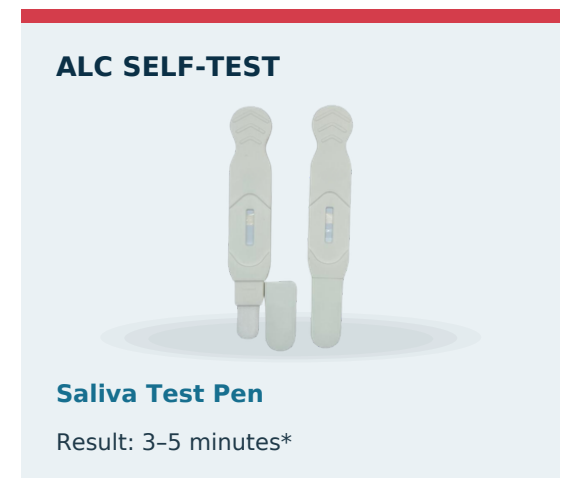

ALC SELF-TEST

Saliva Test Pen

Result: 3-5 minutes*

* According to the catalog. Follow the instructions for the specific test. Photograph is illustrative.

Appendix A | Saliva: Available Substances

33 substances/groups plus alcohol · cutoffs from the supplied list

NR.	SUBSTANCE / GROUP	CUTOFF (ng/mL)
01	Amphetamine (AMP)	50
02	Methamphetamine (MET)	50
03	Cannabis (THC)	50 / 40 / 25 / 20 / 12 / 3,5
04	Phencyclidine (PCP)	10
05	Cocaine (COC)	20
06	Opiates (MOP/OPI)	40 / 25
07	Methadone (MTD)	30
08	MDMA	50
09	Oxycodone (OXY)	50 / 20
10	Cotinine (COT)	20
11	Benzodiazepine (BZO)	50 / 30 / 20 / 10
12	Synthetic Cannabinoids (K2)	25
13	Ketamine (KET)	50
14	Barbiturates (BAR)	50
15	Buprenorphine (BUP)	10 / 5
16	Tramadol (TML)	30
17	6-MAM (Heroin Metabolite)	10

NR.	SUBSTANCE / GROUP	CUTOFF (ng/mL)
18	Fentanyl (FYL)	50 / 20 / 10
19	Carfentanil (CFYL)	50
20	MDPV	300
21	Alpha-PVP	300
22	LSD	10
23	Propoxyphene (PPX)	50
24	Methaqualone (MQL)	300
25	Carisoprodol (CAR)	300
26	EDDP	50
27	AB-PINACA (ABP/K3)	10
28	UR-144 (K4)	25
29	Zopiclone (ZOP)	50
30	Gabapentin (GAB)	2.000
31	Pregabalin (PGB)	500
32	Kratom (KRA)	300
33	Mephedrone (MEP)	100

ALCOHOL (ALC) 0.02 / 0.04 / 0.08 / 0.30% (alcohol; not ng/mL)

Cutoffs vary by product version; confirm availability and the exact cutoff before ordering.

Source: Substante_Dispositive_Saliva_CRP Fusion_02.2026_Social Med_.pdf

Appendix B | Urine: Available Substances

Entries 01–34 · cutoffs from the supplied list

NR.	SUBSTANCE / GROUP	CUTOFF (ng/mL)	NR.	SUBSTANCE / GROUP	CUTOFF (ng/mL)
01	Paracetamol (ACE)	5.000	18	Tramadol (TML)	100 / 200 / 300 / 500
02	Amphetamine (AMP)	1.000 / 500 / 300	19	Ketamine (KET)	1.000 / 500 / 300 / 100
03	Barbiturates (BAR)	300 / 200	20	Oxycodone (OXY)	300 / 100
04	Benzodiazepine (BZO)	500 / 300 / 200 / 100	21	Cotinine (COT)	500 / 300 / 200 / 100 / 50 / 10
05	Buprenorphine (BUP)	10 / 5	22	EDDP	300 / 100
06	Cocaine (COC)	300 / 200 / 150 / 100	23	Fentanyl (FYL)	200 / 100 / 20 / 10
07	Cannabis (THC)	600 / 300 / 200 / 150 / 50 / 30 / 10	24	Synthetic Cannabinoids (K2)	50 / 30 / 25
08	Methadone (MTD)	300 / 200	25	6-MAM	10
09	Methamphetamine (MET)	1.000 / 500 / 300	26	MDA	500
10	MDMA	300 / 500 / 1.000	27	Ethyl Glucuronide (ETG)	1.000 / 500 / 300
11	Morphine/Opiates (MOP/OPI)	300 / 200 / 100	28	Clonazepam (CLO)	400 / 150
12	Methaqualone (MQL)	300	29	LSD	10 / 20 / 50
13	Meperidine (MPRD)	100	30	Methylphenidate (MPD)	300 / 150
14	Opiates (OPI)	2.000 / 1.000	31	Ritalinic Acid (MPD)	1.000
15	Phencyclidine (PCP)	50 / 25	32	Zolpidem (ZOL)	50
16	Propoxyphene (PPX)	300	33	Mephedrone (MEP)	500 / 100
17	Tricyclic Antidepressants (TCA)	1.000 / 500 / 300	34	MDPV	1.000 / 500 / 300

The source lists cutoff variants. Confirm the exact panel and cutoff before ordering.

Source: Substante_Dispositive_Urina_CRP Fusion_02.2026_Social Med_.pdf

Appendix B | Urine: Available Substances

Entries 35–67 · cutoffs from the supplied list

NR.	SUBSTANCE / GROUP	CUTOFF (ng/mL)
35	Diazepam (DIA)	300 / 200
36	Zopiclone (ZOP)	50
37	Methcathinone (MCAT)	500
38	7-Aminoclonazepam (7-ACL)	300 / 200 / 100
39	Carfentanil (CFYL)	500 / 250
40	Cannabinol (CNB)	500
41	Caffeine (CAF)	1.000
42	Cathine (CAT)	150
43	Tropicamide (TRO)	350
44	Alprazolam (ALP)	100
45	Pregabalin (PGB)	50.000 / 500
46	Gabapentin (GAB)	2.000
47	Zaleplon (ZAL)	100
48	Carisoprodol (CAR)	2.000 / 1.000
49	AB-PINACA (ABP)	10
50	Quetiazepam (QTP)	1.000
51	Fluoxetine (FLX)	500

NR.	SUBSTANCE / GROUP	CUTOFF (ng/mL)
52	UR-144	25
53	Kratom (KRA)	300
54	Tilidine (TLD)	50
55	Trazodone (TZD)	200
56	Alpha-PVP	2.000 / 1.000 / 500 / 300
57	Mescaline (MES)	100 / 300
58	Papaverine (PAP)	500
59	Citalopram (CIT)	500
60	Fluoroketamine (FKET)	1.000
61	Olanzapine (OZP)	1.000
62	Risperidone (RPD)	150
63	Tapentadol (TAP)	1.000
64	DMT (NND)	1.000
65	Scopolamine (SCOP)	500
66	Mirtazapine Metabolite (MTZ)	500
67	Alcohol (ALC)	0,02 %

Alcohol is expressed in percent; all other values are in ng/mL. Confirm the version and cutoff.

Source: Substante_Dispozitive_Urina_CRP Fusion_02.2026_Social Med_.pdf

Price List | Saliva, Urine and Drink-Check

2026 price list · EUR net and RON net / gross · source dated 28 Sep 2026

PRODUCT / CODE	EUR	RON	RON
	net	net	gross
ALC Pen with Box CRP-D062-8003-NM	2,66	13,95	16,88
ALC Pen without Box CRP-D062-8003-NM	2,22	11,68	14,13
THC Cassette CRP-D022-P8002-NM	4,33	22,71	27,48
THC Pen CRP-D022-8003-NM	4,66	24,46	29,60
COC Pen CRP-D006-8003-NM	4,66	24,46	29,60
KET Pen CRP-D011-8003-NM	5,99	31,47	38,08
LSD Pen TZKM-D035-8003-NM	6,33	33,22	40,20
PCP Pen TZKM-D019-8002-NM	4,66	24,46	29,60
MDMA Pen CRP-D012-8003-NM	4,33	22,71	27,48
8+1 Cassette CRP-D000-8085-A-NM	7,75	40,71	49,25
Fit2Work 12-Target Pen CRP-D000-8123-D-NM	8,62	45,26	54,77
Production 12-Target Cup CRP-D000-8127PC-NM	8,50	44,61	53,98
Office 12-Target Cup CRP-D000-8127OC-NM	8,60	45,15	54,63
Fit2Drive PRO 15-Target Cup CRP-D000-8157DC-NM	8,09	42,45	51,36

PRODUCT / CODE	EUR	RON	RON
	net	net	gross
5-Target Panel CRP-D000-1054-NM	3,52	18,51	22,39
6-Target Panel CRP-D000-1064-NM	3,91	20,52	24,83
7-Target Panel CRP-D000-1074-NM	4,53	23,78	28,78
8-Target Panel CRP-D000-1084-NM	4,83	25,33	30,65
9-Target Panel CRP-D000-1094-NM	5,15	27,04	32,72
10-Target Panel CRP-D000-1104-NM	5,48	28,75	34,78
12-Target Panel CRP-D000-1124-NM	7,87	41,32	50,00
20-Target Panel CRP-D000-1204-NM	13,40	70,34	85,11
Clinics 10-Target Cup CRP-D000-1107-NM	7,13	37,44	45,30
Wristband 3x3 HD-000-D036-NM	11,74	61,65	74,59

Prices are taken from the supplied list. Confirm price and availability before ordering.

Source: „V2 CRP Online Shop Price List VX 200626.pdf“.

FAQ | Products and Test Formats

WHAT DO THE PRODUCT NAMES MEAN?

Fit2Drive refers to products for personal screening before a planned journey. Fit2Work covers configurations for workplace testing programs. Always check the exact substance list and instructions for the selected product.



TEST FORMATS AT A GLANCE

FORMAT	WHAT IS IT?	EXAMPLE IN THIS BROCHURE
Pen test	A saliva collector with a test for one or more substances.	ALC, THC, Fit2Work
Cassette	Saliva is collected and placed in the test well.	THC 3,5 ng/mL; Fit2Drive 8+1
Cup	The sample is collected and tested in one cup system.	Fit2Drive PRO, 15 substances
Panel	The specific set of substances or substance groups in a test.	Urine panel with 9 substances

8+1 = eight drug groups plus a separate alcohol indicator. Check the product for its exact composition.

FAQ | Results, Cutoffs and Confirmation

What is a screening test?

An initial qualitative check. The result shows whether the test reacts at the criteria set for that product.


What does cutoff mean?

The concentration at which a test is designed to react. A value of 3.5 ng/mL refers to the stated specimen and substance; it is not a blood value.



HOW DO I READ A TEST WITH C AND T LINES?

NEGATIVE

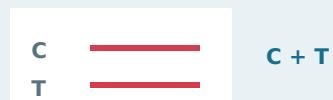

The substance was not detected at the test cutoff.

PRELIMINARY POSITIVE



Seek appropriate confirmation before important decisions.

INVALID



Repeat with a new test, following its instructions.

Detection windows vary. Saliva or urine alone cannot establish when a substance was used or fitness to drive.

Follow the instructions for the specific test. The C/T pattern applies to drug test strips; alcohol indicators may differ.
General reference: FDA, „Drugs of Abuse Home Use Test” ([fda.gov/medical-devices/drugs-abuse-tests/drugs-abuse-home-use-test](https://www.fda.gov/medical-devices/drugs-abuse-tests/drugs-abuse-home-use-test)).

FAQ | PEDT Services and Test Selection



What is PEDT?

Pre-Employment Drug Testing: drug testing arranged before hiring, under agreed requirements and applicable procedures.

What is background screening?

Verification of information relevant to hiring. This separate service can be coordinated with PEDT.

PEDT SERVICE WORKFLOW

01 Request

Country, location, panel



02 Appointment

Test center and candidate



03 Sample collection

Sample and identity



04 Screening

Preliminary result



05 Report

As agreed with the client

What is random testing?

Testing under a predefined company policy, using specified intervals or selection procedures.

What should an order specify?

Specimen: saliva or urine; substances, cutoffs, quantity and location. Request the data sheet for the exact variant.